


MRF-25V300

Combined Refrigerator & Freezer

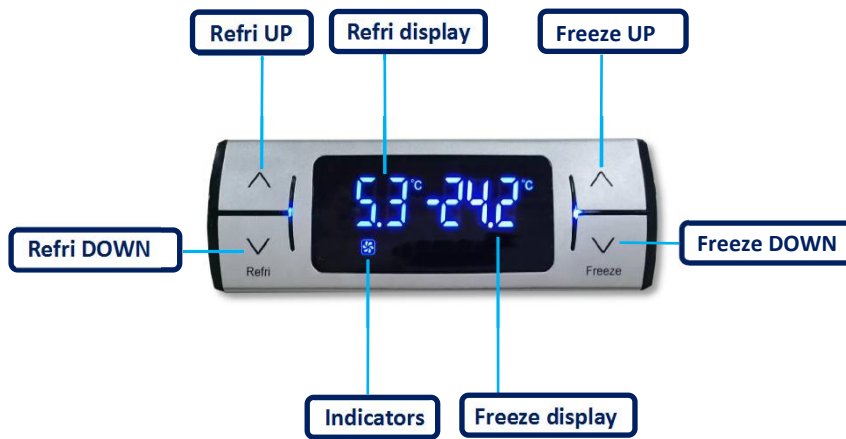

MRF-25V300 Medical Combined Refrigerator and Freezer (upper 2~8°C and lower -25°C) is widely used in hospitals, blood stations and laboratories. The product has double chamber, double display, double control system, two sets of individual compressor and refrigeration system. Temperature in both chambers can be independently adjusted.

FEATURES



Multi-Scene applications
Double chamber & smaller footprint



Heated glass design
Fast defogging & no condensation



Coil evaporator
Better performance, easy maintenance



Suction circulation system
Less temp. fluctuation & better



Lockable doors
Standard locks & padlocks



Large LED displays
Easy to observe and operate, built-in USB



Precise temp. control
Temp. control accurate to 0.1°C



6 types of alarm
Audible and visual alarms



Variable space
Adjustable shelves and drawers




Eco-friendly
CFC-free refrigerant and energy-saving

STRUCTURE

Optional Accessories

Data logger
RS232 interface
RS485 interface

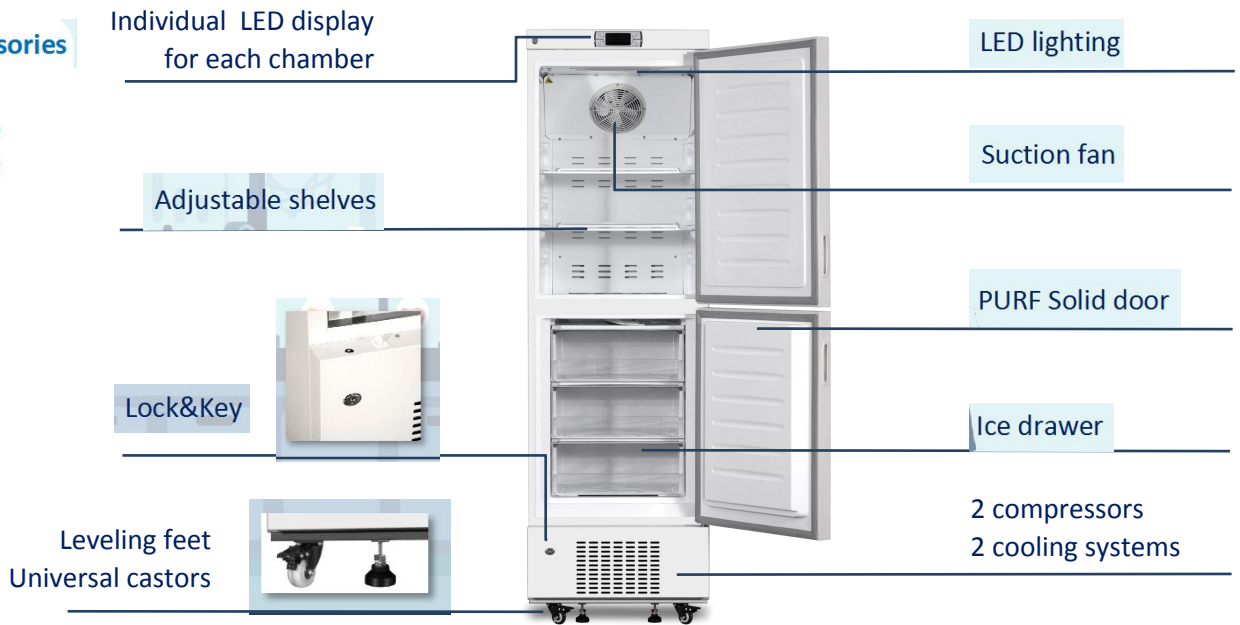


Combined R&F

High/low temp. alarm

Power failure alarm

Door ajar alarm

Thermostat failure alarm

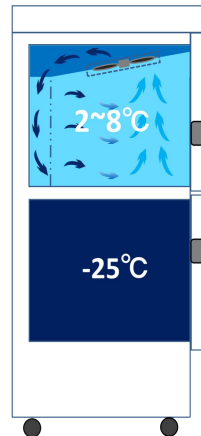

Refrigerator air-force circulation system keeps the chamber temperature fluctuation extremely stable at $5 \pm 3^{\circ}\text{C}$.

The coil evaporator evenly distributed around the cabinet. Perfect control of the temperature uniformity in the chamber

Side cross-sectional view

PARAMETER

Performance		Refrigeration System		Alarms	
Cooling method	Direct cooling	Capacity (L/cu.ft)	R:198 F:102	High/low temperature	Y
Refrigerator Defrost	Auto	Temp. controller	Microprocessor	Power failure	Y
Freezer defrost	Manual	Sensor	NTC	Sensor error	Y
Temperature range	F:-10°C~-25°C R:2°C~8°C	Display	Digital display	Controller failure	Y
Ambient temperature	10~32°C	External material	SPCC sprayed steel	Door ajar	Y
Noise (db)	49	Internal material	HIPS	Ambient temperature	/
Electrical Data		Dimensions		Other	
Voltage frequency(V/Hz)	220V 50Hz/60Hz	Package size(W*D*H)(mm)	720*695*2015	Inside door	3 shelves & 6 drawers
Power (W)	349W	NT./GT. (kg)	100/120	Remote alarm port	Optional
Current (A)	2.2A	Exterior size(W*D*H)(mm)	650*620*1980	Access port QTY/Dia.	2/25MM
Power consumption (kWh/24h)	2.96	Interior size(W*D*H)(mm)	R:565*480*759 F:461*465*600	Data logger	Optional