

# SafeVac Vacuum-Controlled Aspiration System

Safe and efficient system for labs handling critical or biological liquid waste.

## Features

- 1.Compact, user-friendly design
- 2.Level detection system prevents bottle overfill
- 3.Control knob for vacuum setting
- 4.The parts in contact with liquids are autoclavable
- 5.Easy to disassemble and clean
- 6.Ideal tool to avoid lab air contamination and protect personnel safety
- 7.Adjustable vacuum level with indicator light is designed for various

## Applications

For laboratory waste recovery and separation of liquid and solid.

It can be widely used in cell culture, DNA extraction, microplate waste removal, and any other liquid separation or recovery.


Best replacement for the common vacuum pump, for all kinds of waste liquid recovery, especially biological, chemical and radioactive liquids.

## Specifications

Wide range of adapters to fit well with different tubes, microplate and disk plate, all of which can be fully autoclaved (At 121°C 1 bar for 20 minutes. After autoclaving the equipment must be cooled down before use)

<b>Vacuum range</b>	<b>0-600mbar</b>
<b>Aspirating speed</b>	<b>15L/min (air)</b>
<b>Flow rate [liquid]</b>	<b>17mL/s</b>
<b>Bottle</b>	<b>4L</b>
<b>Autoclave component</b>	<b>Tube;waste bottle;operation handle;</b>
<b>Noise level</b>	<b>≤55dB</b>
<b>Input voltage</b>	<b>100-240VAC, 50/60Hz</b>
<b>Dimension[WxDxH]</b>	<b>180x240x450mm</b>
<b>Weight</b>	<b>3kg</b>

<b>A</b>	<b>8-channel tip detrusion adapter 200μL tip</b>
<b>B</b>	<b>8 channel needle Ø1.5mm</b>
<b>C</b>	<b>Single needle Ø2.5mm, 40mm</b>
<b>D</b>	<b>Single needle Ø2.5mm, 120mm</b>
<b>E</b>	<b>Single needle Ø1.5mm, 40mm</b>
<b>F</b>	<b>Single needle Ø1.5mm, 120mm</b>
<b>G</b>	<b>tip detrusion adapter 200μL tip</b>
<b>H</b>	<b>tip detrusion adapter 1mL tip</b>
<b>I</b>	<b>SafeVac Handle bracket</b>

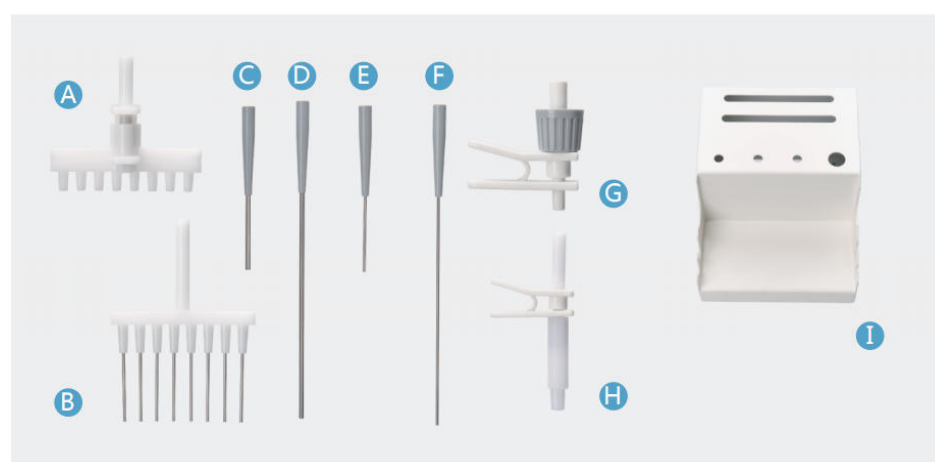

Self-closing connectors to avoid escape of drops or aerosols from disconnected bottle

Stand for transferring the bottle securely and holding the hand operator

vacuum level indicator light

