





MDF-40V528 Medical Deep Freezer (-20°C ~ -40°C) is widely used in hospitals, blood stations and laboratories. The product has double chamber, double display, double control system, two sets of individual compressor and refrigeration system. Both chambers can be independently

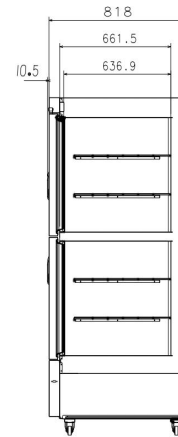

FEATURES

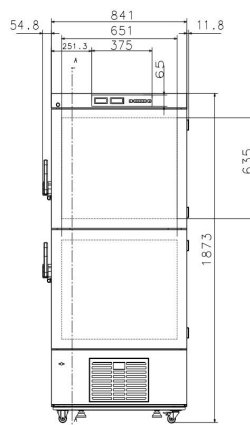

Imported brand compressor
Double cooling system, double chamber



Large LED displays
Easy to observe and operate, built-in U



Fast pull down
150 minutes from +25°C to -40°C



Precise temp. control
Temp. control accurate to 0.1°C



Coil evaporator
Better performance, easy maintenance



6 types of alarm
Audible and visual alarms



Thick PURF insulation
Better insulation performance




Variable space
Adjustable shelves and drawers



Lockable doors
Standard locks & padlocks





Eco-friendly
CFC-free refrigerant and energy-saving

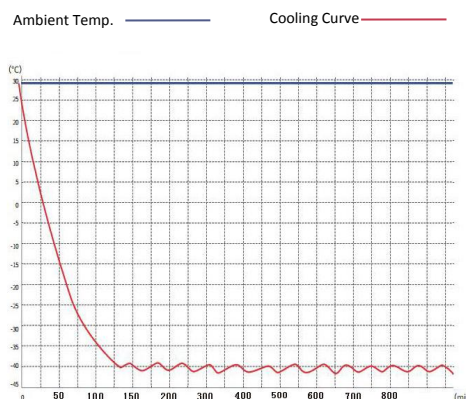

MDF-40V528 Pull down at ambient temp. 32°C

STRUCTURE

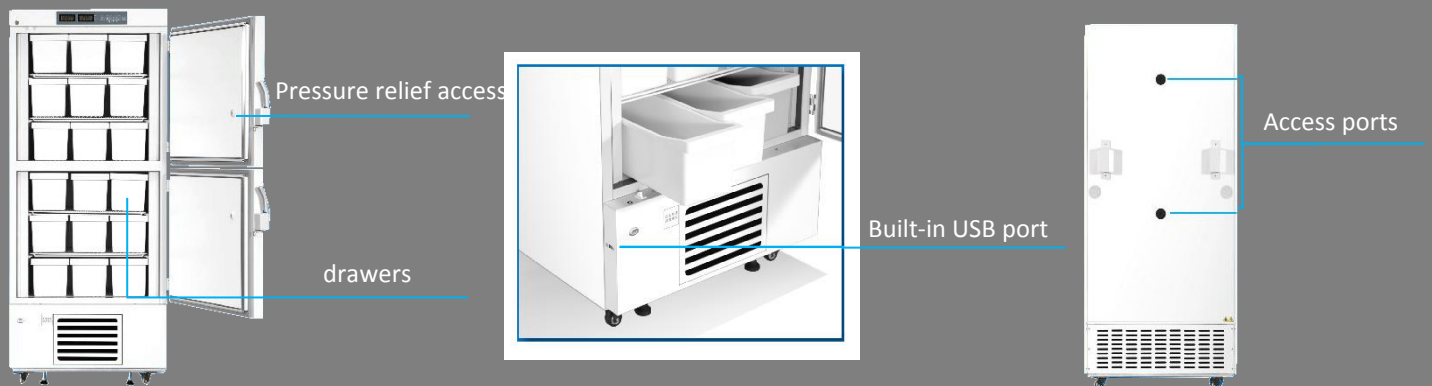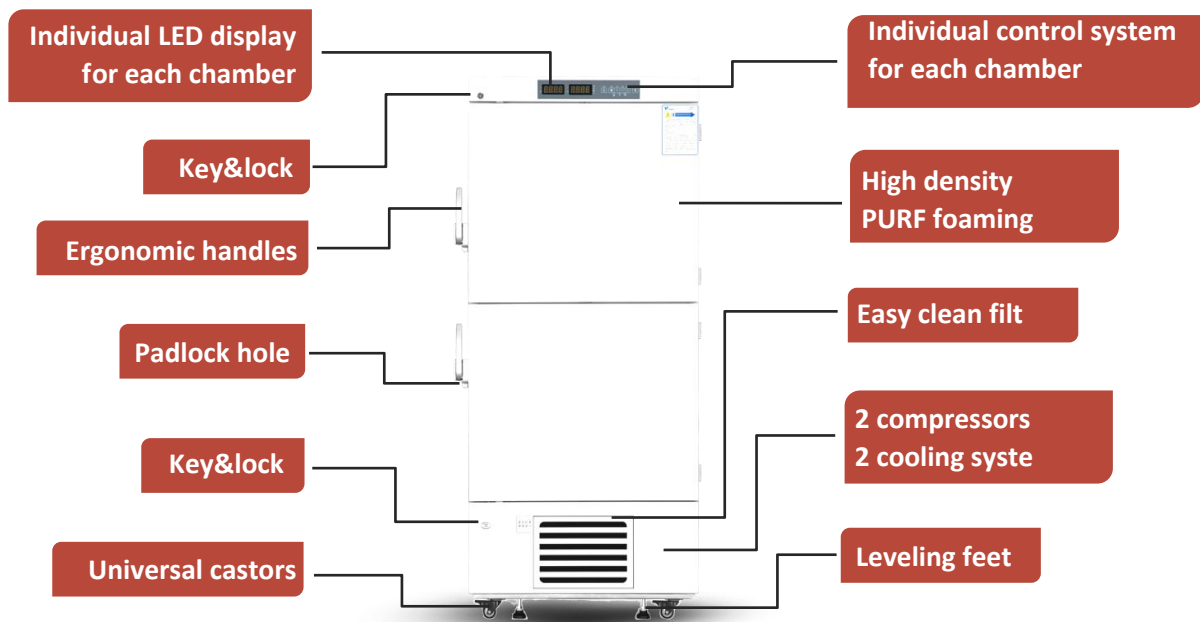

PARAMETER

Performance		Refrigeration System		Alarms	
Cooling method	Direct cooling	Compressor brand/QTY	CUBIGEL/2	High/low temperature	Y
Defrost	Manual	Condenser	Copper tube	Power failure	Y
Refrigerant	R290	Temp.controller	Microprocessor	Sensor error	Y
Temperature range	-20℃~-40℃	Sensor	NTC	Controller failure	Y
Ambient temperature	10~32℃	Display	Digital display	Door ajar	Y

Electrical Data		Dimensions		Other	
Voltage frequency	220/50Hz	Capacity(L/cu.ft)	528/18.64	Remote alarm port	Y
Power(W)	340W	Exterior	910*845*1875	Remote alarm port	Y
Current(A)	3.4 A	Interior	(651*631.5*635)*2	Access port QTY/Dia.	2/25MM