

MIT300/500 Metallurgical Microscope

- Dark field objective, Make observation in dark & bright field.
- Oversize view field eyepiece, view field up to 22mm for comfortable observation.
- Trinocular head with dual transform.
- Light distribution (both): 100: 0 (100% for eyepiece) ; 80:20 (80% for trinocular head and 20% for eyepiece)



MIT300(Led lamp version)



MIT500(Led lamp version)



MIT300(halogen lamp version)




MIT500(halogen lamp version)



- Stage with fast lowering unit can increase 50mm for working distance, more convenient for observing big size sample.
- Oversized three-layer mechanical stage provides more choice for samples.

- EPI reflection illumination device(LED lamp version)

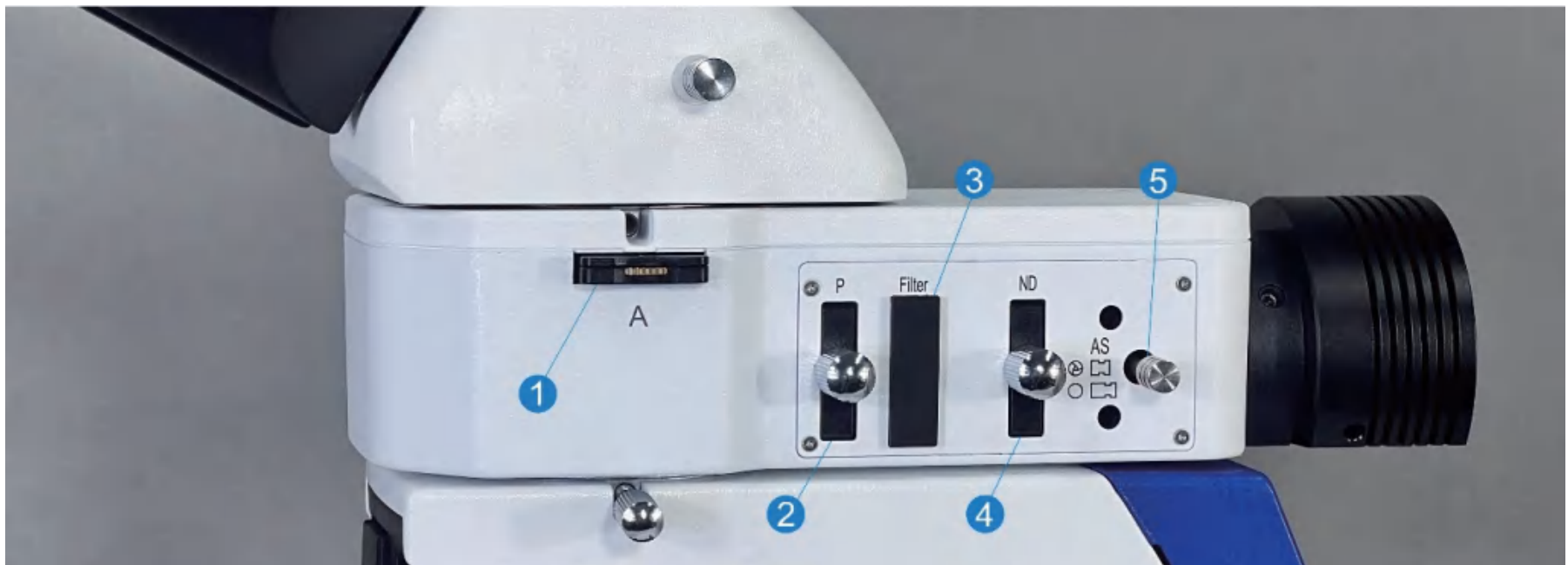

1. Analyzer 2. Polarizer 3. Filter switch handle 4. Ground glass slider 5. Aperture diaphragm

- EPI reflection illumination device(Halogen lamp version)



1. Polarizer 2. Filter switch handle 3. Field diaphragm 4. Aperture diaphragm A. Analyzer B. Ground glass slider

● Objective Specification: ("●" is standard outfits, "○"is optional accessories.)
 “T” means transmitting “R” means reflecting “T&R “ means transmitting and reflecting

Item	Objective	N.A.	W.D.(mm)	MIT300/MIT500	MIT300D/MIT500D
Metallurgical LWD Infinity Plan Objectives for Bright/Dark Field	M PLAN 5× BD	0.15	7.30		●
	M PLAN 10× BD	0.25	11.50		●
	M PLAN 20× BD	0.40	5.60		●
	M PLAN 50× BD	0.75	1.60		●
Metallurgical LWD Infinity Plan Objectives	LMPlan FL 5× (Semi Apochromatic Objective)	0.15	21.00	●	
	LMPlan FL 10× (Semi Apochromatic Objective)	0.30	20.00	●	
	LMPlan FL 20× (Semi Apochromatic Objective)	0.40	15.00	●	
	LMPlan FL 50× (Semi Apochromatic Objective)	0.75	4.25	●	
	LPL 80×	0.80	0.85	○	
	LMPlan 100× (Dry) (Semi Apochromatic Objective)	0.85	3.00	○	




Metallurgical LWD Infinity Plan Objectives




Metallurgical LWD Infinity Plan Objectives for Bright & Dark Field

MIT300 Metallurgical Microscope Size(mm)

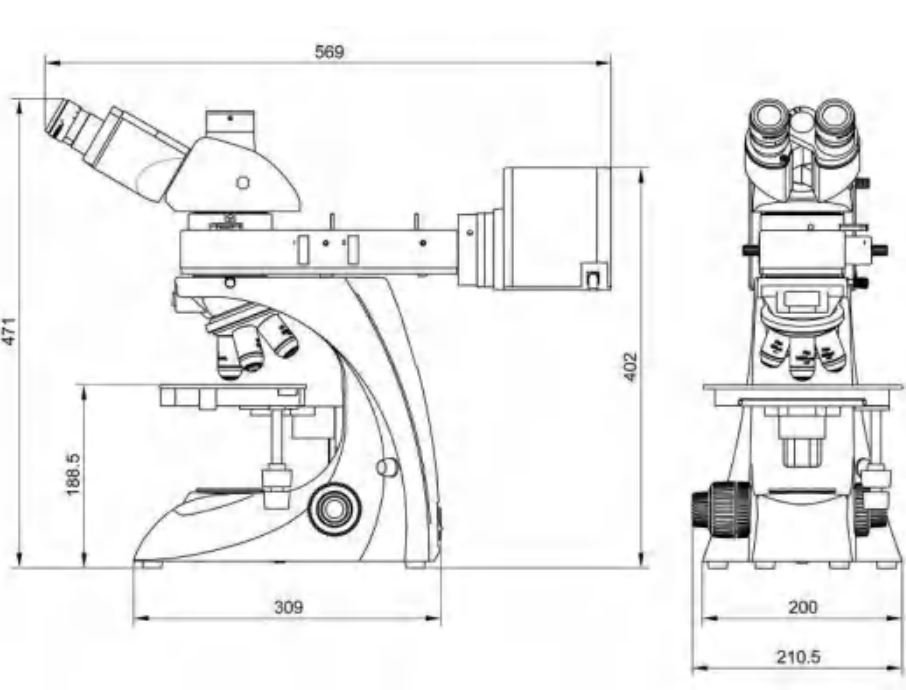


Haolgen Version

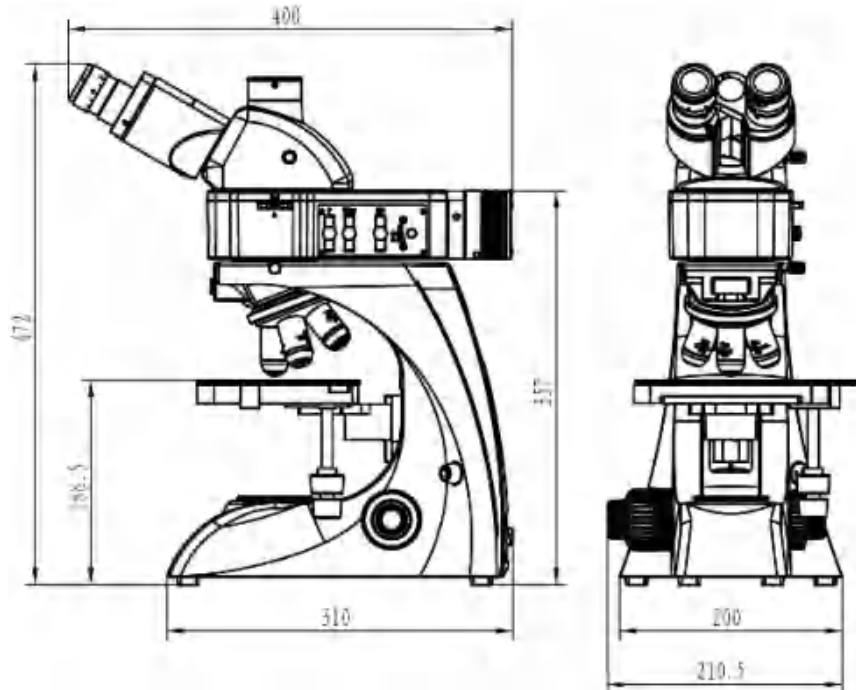


LED Version


MIT500 Metallurgical Microscope Size(mm)

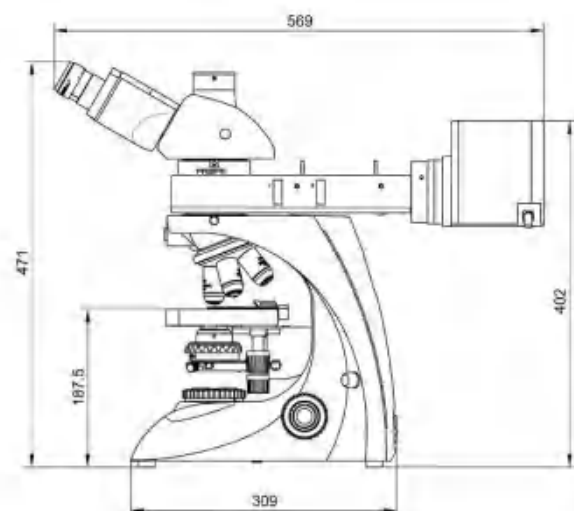


Halogen version

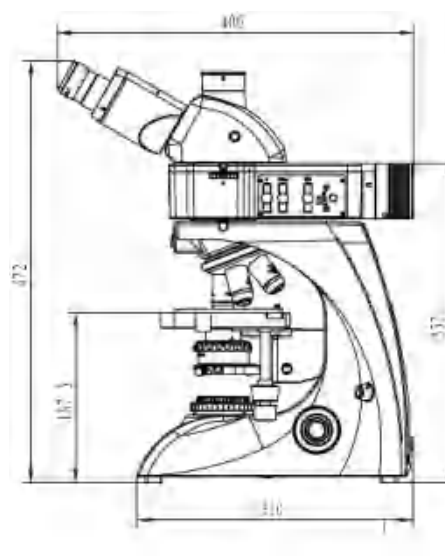

LED version

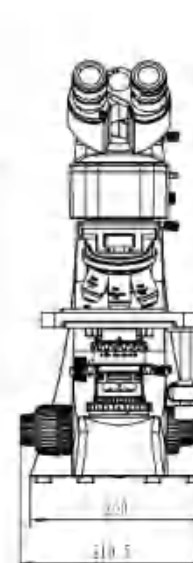

MIT300/500 Metallurgical Microscope General Layout Diagram

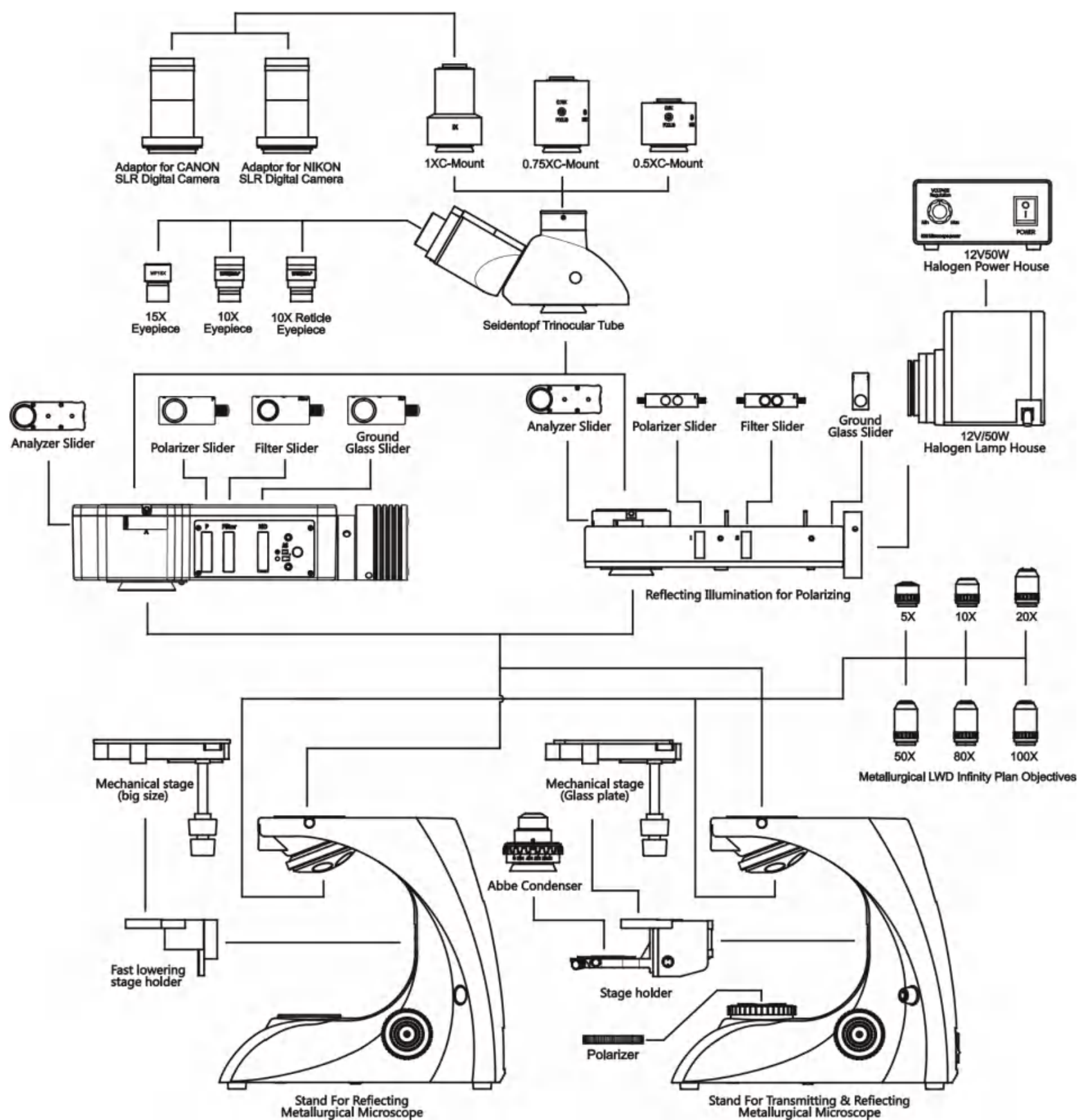

MIT300/500 Metallurgical Microscope Outfits

Item	Specification	MIT300	MIT300D	MIT500	MIT500D
Eyepiece	WF10X／22mm (adjustable)	●	●	●	●
	WF10X／22mm (adjustable , reticule 0.1mm)	●	●	●	●
Metallurgical LWD Infinity Plan Objective	LMPlan FL 5×/0.15 (Semi Apochromatic Objective) W.D=21.00mm	●		●	
	LMPlan FL 10×/0.30 (Semi Apochromatic Objective) W.D=20.00mm	●		●	
	LMPlan FL 20×/0.40 (Semi Apochromatic Objective) W.D=15.00mm	●		●	
	LMPlan FL 50×/0.75 (Semi Apochromatic Objective) W.D=4.25mm	●		●	
	LPL80X／0.80 W.D=0.85mm	○		○	
	LMPlan FL 100×/0.85 (Dry) (Semi Apochromatic Objective) W.D=3.00mm	○		○	
Metallurgical LWD Infinity Plan Objectives for Bright & Dark Field	M PLAN 5×／0.15 BD W.D=7.30mm		●		●
	M PLAN 10×／0.25 BD W.D.= 11.50mm		●		●
	M PLAN 20×／0.40 BD W.D.= 5.60mm		●		●
	M PLAN 50×／0.75 BD W.D.= 1.60mm		●		●
Seidentopf Trinocular Head	Inclined 30°,rotatable 360°, Interpupillary distance: 48mm～76mm Light distribution (both): 100: 0 (100% for eyepiece) 80:20 (80% for trinocular head and 20% for eyepiece)	●	●	●	●
Nosepiece	Quintuple	●		●	
	Quadplex (specially for bright /dark field objectives M26x1/36")		●		●
Mechanical Stage	Stage size: 200mm x 150mm Travel range: 77mm x 52mm Stage fast lowering unit: Up-Down moving range can be increased 50mm.	●	●		
	Stage size: 182mm x 140mm, Travel range: 77mm x 52mm			●	●
Polarizing Unit	Polarizer /Analyzer	●	●	●	●
Condenser	Abbe N.A.1.25 , with Iris Diaphragm.			●	●
Coaxial Coarse and Fine Focusing Unit	Coarse stroke: 25mm, Fine stroke per rotation: 0.2mm, fine division 2 μm, Tightness adjustment, upper distance limiter.	●	●	●	●
Reflecting Illumination	5WLED lamp DC input (Universal voltage: 100V～240V)	●		●	
	12V/50W Halogen lamp (Universal voltage: 100V～240V)	○	●	○	●
	10WLED lamp DC input (Universal voltage: 100V～240V)			○	○
Transmitting Illumination	5WLED lamp DC input (Universal voltage: 100V～240V)			●	
	12V/30W Halogen lamp (Universal voltage: 100V～240V)			○	●
	10WLED lamp DC input (Universal voltage: 100V～240V)				○
Filter for transmission lamp(Φ45)	Blue			●	●
	Green/ Amber/ Grey			○	○
Filter for Reflection lamp(Φ21)	Blue	●	●	●	●
	Green/ Amber/ Grey	○	○	○	○
C-Mount	1 X /0.75X/0.5X/0.35X(Focus Adjustable)	○	○	○	○

Note: "●"In Table Is Standard outfits, "○" Is Optional Accessories.

Packing Size: 520mm×470mm×450mm Gross Weight: 15.5 kgs Net Weight: 14kgs