

## DW-LAC-N-II

# Climate Chamber with Vertical Light

### Features

- ★ Two sets of refrigeration systems automatically switch in turn to ensure that long time operation without failure.
- ★ International brand fluorine-free compressor, environment friendly and reliable.
- ★ International brand humidity sensor, high accuracy for humidity.
- ★ Double door design. The inner door is made of tempered glass for easy observation, the outer door with magnetic strip, good sealing.
- ★ Large LCD screen to display temperature and humidity and light intensity. Multi-segment timing program control, high precision.
- ★ Polished stainless-steel chamber, moveable grids can be freely adjusted, easy for cleaning. Powder coating steel housing.
- ★ External water tank. Electric heater humidification system.
- ★ Two adjustable LED light boards, LED stepless dimming.
- ★ The parameters can be saved automatically in case of power failure, and the original setting programme will be run when power is restored.
- ★ In case of failure of the incubator, a fault message appears on the LCD display, and the operational faults can be seen at a glance.
- ★ JAKEL circulation fan, unique air duct structure, effectively guarantee the temperature stability inside the box.
- ★ Equipped with leakage protection and a spare mechanical thermostat to ensure the product work normally even the main temperature control failed (for heating).
- ★ USB port for recording data.


### Safety Device

1. Over temperature protection
2. Compressor over current protection
3. Over compressing protection
4. Water shortage protection
5. Heater overheats protection
6. Audible and visual alarm system

### Options

CO2 control system

White/Blue/Red/Yellow/Plant growth light board

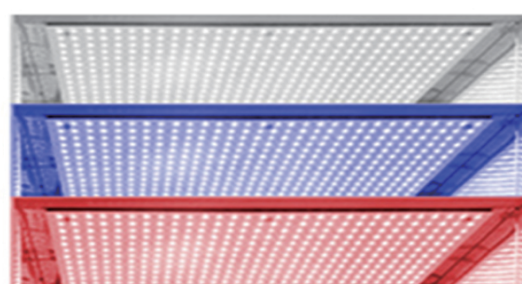

The light plate has red light, blue light, white light, plant growth light to choose from.



## Specifications

Model	DW-LAC-175-N-II	DW-LAC-275-N-II	DW-LAC-375-N-II	DW-LAC-475-N-II	DW-LAC-800-N-II	DW-LAC-800D-N-II	DW-LAC-1075-N-II	DW-LAC-1075D-N-II	DW-LAC-1675D-N-II
Chamber Volume	175L	275L	375L	475L	800L	800L	1075L	1075L	1675L
Temp. Range	With light:10~55℃ Without light:0~55℃								
Temp. Stability	≅±1℃								
Display Resolution	0.1℃								
Temp. Uniformity	≅±1℃								
Humidity Range	30~95%RH								
Humidity Deviation	±3%RH								
Light Intensity (LX)	2 LED boards White light: 0-20000LX stepless dimming Other lights: 0-5000LX stepless dimming								
Power Supply	AC220V, 50HZ								
Operating Temp.	+0~40℃								
Chamber Size(W×D×H)cm	45*42*93	58*51*94	59*55*116	70*55*125	96*61*137	126*50*127.5	95*70*160	126*67*127.5	146*85*135
Exterior Size(W×D×H)cm	65*73*175	77*74*171	78*87*191.5	88*87*199.5	110*93*227	141.8*78.1*195	101*90*234	141.8*95.1*195	162.5*122*204
Grids(Std/Max)	3/16	3/16	3/21	3/22	3/22	6/18(2 grids in one level)	3/22	6/18(2 grids in one level)	6/20(2 grids in one level)

D in the model# means double doors.